


A1c Go

Boronate Affinity Chromatography POCT HbA1c Analyzer




Compact

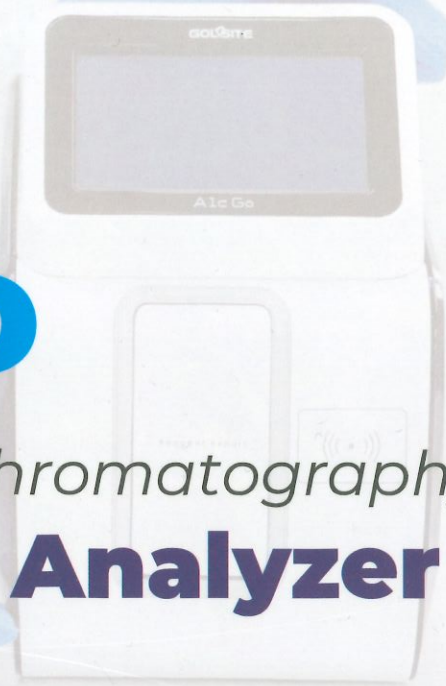

Fully Automatic



Cartridge based

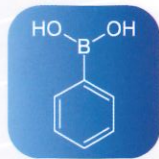


Boronate Affinity



Features


- Advanced technique: Boronate Affinity Chromatography
- Fully automatic one-step operation
- Fast and accurately answer of results in minutes
- Least interference from HbA1c Variants
- No liquid Path: no carryover & minimized maintenance
- Point-of-Care analyzer is suitable for deliver timely primary care at the physician's office

Smart Look

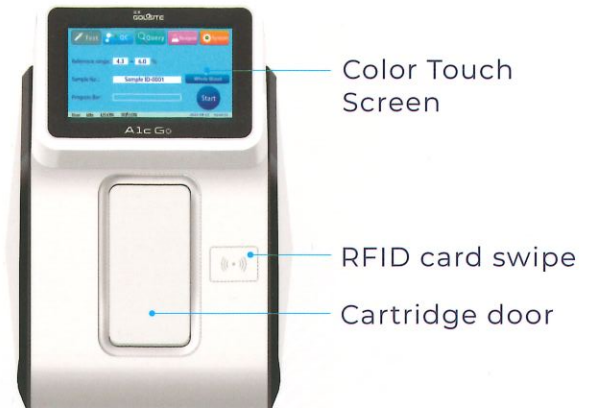

224mm(L)*160mm(W)*236mm(H)


Simplified Operation

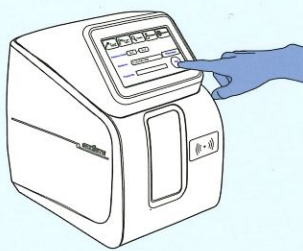

Click to START



Collect 5 μ L sample, insert the pipette into R1 position

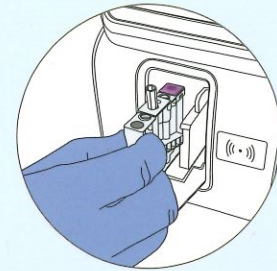

Insert the cartridge into the machine, followed by automatically testing process. Result will be displayed in 5 min

Specifications

Screen	Color touch screen	Methodology	Boronate Affinity Chromatography
Assay	HbA1c, Glycated Hemoglobin	Pack size	25 kits/box
Sample volume	5 μ L	Linear range	3.5-15.0% (NGSP)
Interval	\leq 5 min	Calibration	RFID card swipe
CV	\leq 3%	Dimensions	224mm(L)*160mm(W)*236mm(H) \leq 5 Kg
Communications	LAN port, USB, Serial Port, WiFi		
consumable	Cartridge, RFID card, disposable pipette, HbA1c Control(optional)		

Goldsite Diagnostics Inc.

ADDRESS: No. 103C, 503C & 504D, Technology Building & No. 3A & 4A, Technology Building Annex, Zhaoshang Sub-District, Nanshan District, Shenzhen, China, 518067
 TEL: 0086 755 2689 0807 FAX: 0086 755 2689 0799 EMAIL: sales@goldsite.com.cn
 WEBSITE: www.goldsite.com.cn



@ Goldsite Diagnostics